

AM7515MT2A



Brand:
Product Code: AM7515MT2A

Short Description

5 Megapixel resolution
130-220 x magnification
Coaxial and brightfield illumination
Flexible LED Control (FLC)
Automatic Magnification Reading (AMR)

Description

The AM7515MT2A is a 5MP Dino-Lite Edge Series microscope that has a magnification range of 130-220x with built-in coaxial illumination. By using the Flexible LED Control (FLC), this model is capable of switching between or mixing, brightfield and coaxial illumination.



USB 2.0



Measurement
functionality



Adjustable
~130x-220x




Metal
housing



Flexible LED
Control FLC



Automatic
Magnification
Reading (AMR)

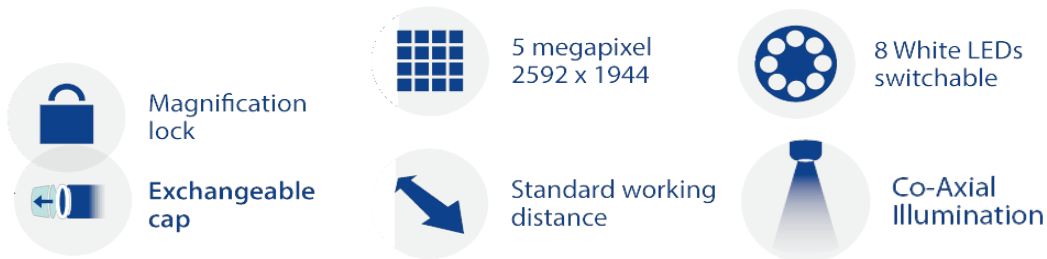

The coaxial illumination technique reveals details that are very difficult to see under normal light. The AM7515MT2A includes the Automatic Magnification Reading (AMR) feature which automatically detects and displays the magnification in the included Dino-Lite software.

These unique features make the Dino-Lite AM7515MT2A a great inspection tool for material analysis, electronics inspection, or any similar application that requires coaxial illumination, versatility and mobility. The Dino-Lite AM7515MT2A is bundled with the user-friendly DinoCapture 2.0 software for Windows.

For this model it includes functions such as:

- Automatic Magnification Reading (AMR)
- Flexible LED Control (FLC)
- Calibration
- Measurement
- Capturing & annotating images
- And recording video

Working distance/field of view/depth of field

| MAGNIFICATION RATE | WORKING DISTANCE* | FIELD OF VIEW(X) | FIELD OF VIEW(Y) | DEPTH OF FIELD |
|---------------------------|--------------------------|-------------------------|-------------------------|-----------------------|
| 130 | 4.5 | 3.0 | 2.2 | |
| 140 | 5.0 | 2.8 | 2.1 | |
| 150 | 5.6 | 2.6 | 1.9 | |
| 160 | 6.3 | 2.4 | 1.8 | |
| 170 | 7.1 | 2.3 | 1.7 | |
| 180 | 8.0 | 2.2 | 1.6 | |
| 190 | 8.9 | 2.1 | 1.5 | |
| 200 | 9.9 | 2.0 | 1.5 | |

| | | | | |
|--|---------------------------|-----|-----|------------------|
| 210 | 10.9 | 1.9 | 1.4 | |
| 220 | 11.9 | 1.8 | 1.3 | 0.1 |
| Listed values may differ slightly | *Without front cap | | | Unit = mm |

Specification

| Lighting | |
|-----------------------------------|---|
| Light/ LED type | White |
| Number of LEDs | 8 |
| LED on/off switchable: | Yes |
| Infrared filter | IR cut-filter >650 nm |
| Diffuser available | Yes, N3C-D |
| Emission filter | No |
| Polarizer | No |
| Optics | |
| Magnification | 130-220x |
| Macro zoom | No |
| Working distance | Standard |
| Lens type | Glass with anti-reflection coating |
| Sensor | |
| Sensor type | CMOS |
| Resolution | 5 Megapixel (2592x1944) |
| Maximum frame rate | 30fps (max 20fps video recording) |
| Compatibility | |
| Interface | USB 2.0 |
| Operating system | Windows 7, 8, 10 & 11, MacOS 10.12 and up |
| Software | DinoCapture 2.0 (Windows), DinoXcope (Mac OS) |
| Supported image formats (Windows) | BMP, GIF, PNG, JPG, TIF, RAS, PNM, TGA, PCX, MNG, WBMP, JP2, JPC, PGX |
| Supported video formats (Windows) | WMV, FLV, SWF |
| Supported image formats (MacOS) | JPEG, PNG |
| Supported video formats (MacOS) | MOV (max 1.3MP) |
| Imaging standards | DirectShow, UVC |
| Wifi | Wireless-ready, requires the WF-10 WiFi streamer (optional) |
| Housing | |
| Housing material | Metal housing |

| | |
|----------------------|---|
| Magnification lock | Yes |
| Dimensions | 10.7cm (L) x 3.2cm (D) |
| Weight | 139g |
| Cable length | 1.8m |
| Features | |
| Special feature | Automatic Magnification Reading (AMR), Flexible LED Control (FLC), Coaxial illumination |
| Measurement | Yes |
| Calibration | Yes |
| Microtouch sensor | Yes |
| ESD safe | Yes |
| Information | |
| Package contents | Microscope, Carry pouch, Software CD, User manual, Quick guide, Calibration target, Front cover box |
| Warranty information | 2 years European warranty |
| Regulatory approval | CE, FCC, ROHS |
| Price range | €1000,00 - €1200,00 |

Product Gallery