

# AM4013MTL-FVW

**Brand:**

**Product Code:** AM4013MTL-FVW



## Short Description

1.3 Megapixel resolution  
Robust metal housing  
395nm + white LEDs  
10-90x magnification  
Long working distance

## Description

The VW series has two types of LED lights that are switchable between white and UV light. This model has a long working distance and a filter to filter out the returning UV light so that an image with natural colors can be obtained. The aluminium housing offers enhanced protection and durability and has the best look and feel.



USB 2.0



Measurement  
functionality



Continuous  
~10-90x



1,3 Megapixels  
1280 x 1024



UV-LEDs  
395nm



Switchable  
UV/white LEDs



Metal housing




Long working distance

## Working distance/field of view/depth of field

MAGNIFICATION RATE	WORKING DISTANCE	FIELD OF VIEW(X)	FIELD OF VIEW(Y)	DEPTH OF FIELD
<b>Listed values may differ slightly</b>				<b>Unit = mm</b>
20	151.0	19.6	15.6	6
30	101.0	13.0	10.4	3.4
40	72.2	9.8	7.8	2.2
50	64.1	7.8	6.3	1.65
60	56.2	6.5	5.2	-
70	51.3	5.6	4.5	-
80	48.3	4.9	3.9	-
90	46.6	4.3	3.5	0.8

## Special Light

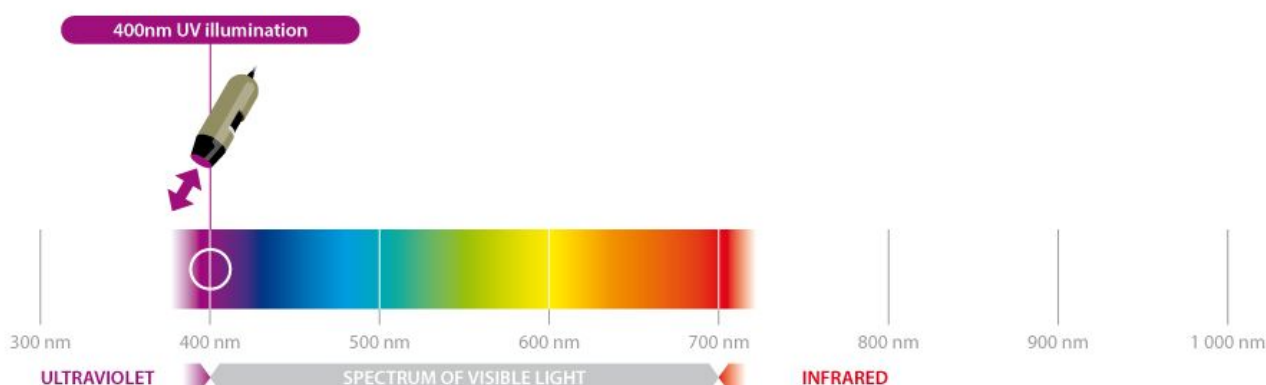

## Specification

Lighting	
Light/ LED type	395 nm UV + white
Number of LEDs	4+4
LED on/off switchable:	Yes
Infrared filter	IR cut-filter >650 nm
Diffuser available	No
Emission filter	Yes, UV cut filter
Polarizer	No
Optics	

Magnification	10x ~ 90x
Macro zoom	No
Working distance	Long
Lens type	Glass with anti-reflection coating
Sensor	
Sensor type	CMOS
Resolution	1.3 Megapixel (1280x1024)
Maximum frame rate	30 fps
Compatibility	
Interface	USB 2.0
Operating system	Windows 7, 8, 10 & 11, MacOS 10.9 and up
Software	DinoCapture 2.0 (Windows), DinoXcope (Mac OS)
Supported image formats (Windows)	BMP, GIF, PNG, JPG, TIF, RAS, PNM, TGA, PCX, MNG, WBMP, JP2, JPC, PGX
Supported video formats (Windows)	WMV, FLV, SWF
Supported image formats (MacOS)	JPEG, PNG
Supported video formats (MacOS)	MOV
Imaging standards	DirectShow, UVC
Wifi	Wireless-ready, requires the WF-10 WiFi streamer (optional)
Housing	
Housing material	Metal housing
Magnification lock	No
Dimensions	10.5cm (L) x 3.2cm (D)
Weight	140g
Cable length	1.8m
Features	
Special feature	No
Measurement	Yes
Calibration	Yes
Microtouch sensor	Yes
ESD safe	Yes
Information	
Package contents	Microscope, carry pouch, software CD, calibration target, user manual
Warranty information	2 years European warranty
Regulatory approval	CE, FCC, ROHS
Price range	€350,00 - €550,00

**Product Gallery**