

AM7915MZTL - EDGE

Brand:
Product Code: AM7915MZTL



Short Description

- 5 Megapixel resolution
- Extended Depth of Field (EDOF)
- Extended Dynamic Range (EDR)
- Flexible LED Control (FLC)
- 10-140x magnification
- Long working distance (LWD)

Description

The AM7915MZTL is part of the Long Working Distance (LWD) range and with its Long Working Distance and magnification range from 10 - 140x, it is suitable for a wide range of applications where extra distance to the object or a larger field of view is required.



USB 2.0



Continuous
~10x - 140x



5 megapixel
2592 x 1944




Long working
distance



Metal
housing



Extended
Depth Of Field
EDOF

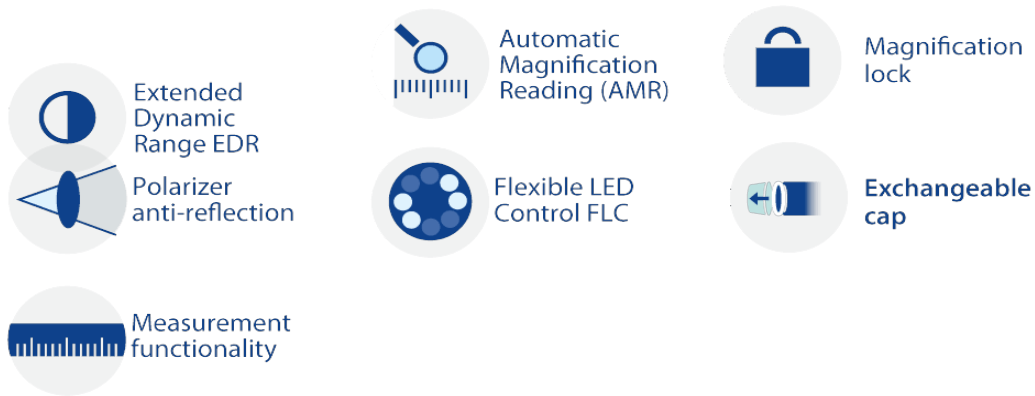

With the Extended Dynamic Range (EDR) feature, the details of darker or brighter areas within the object can be revealed by stacking images at different exposure levels. The Extended Depth of Field (EDOF) feature automatically stacks images at different focus level to improve the depth of field on rough or uneven surfaces. The AMR function is of particular importance when the measurement functions are frequently used. Because of the built-in polarization filter this model is ideal when working with shiny or reflective objects such as metal, plastic, glass, jewelry, electronics, etc.

The main features of the AM7915MZTL are:

- 5 Megapixel Edge sensor
- Extended Depth of Field (EDOF)
- Extended Dynamic Range (EDR)
- 10-140x magnification
- Long working distance
- Automatic Magnification Reading (AMR)
- Adjustable polarizer
- Flexible LED Control (FLC)
- And more...

*(*EDOF/EDR only functions under Windows OS)*

Working distance/field of view/depth of field

MAGNIFICATION RATE	WORKING DISTANCE*	FIELD OF VIEW(X)	FIELD OF VIEW(Y)	DEPTH OF FIELD
Listed values may differ	*Without front cap			Unit = mm

slightly

10	234.5	39.0	29.2	12.0
20	112.5	19.5	14.6	4.5
30	72.5	13.0	9.7	3.1
40	54.5	9.8	7.3	1.8
50	43.5	7.8	5.8	1.2
60	37.5	6.5	4.8	0.9
70	33.7	5.6	4.2	0.6
80	31.3	4.9	3.6	0.5
90	29.9	4.3	3.2	0.37
100	29.2	3.9	2.9	0.28
110	29.1	3.5	2.6	0.22
120	29.3	3.3	2.4	0.17
130	29.8	3.0	2.2	0.12
140	30.5	2.8	2.1	0.09

Specification

Lighting	
Light/ LED type	White
Number of LEDs	8
LED on/off switchable:	Yes
Infrared filter	IR cut-filter >650 nm
Diffuser available	Yes (N3C-D included)
Emission filter	No
Polarizer	Yes, linear
Optics	
Magnification	10x ~ 140x
Macro zoom	No
Working distance	Long
Lens type	Glass with anti-reflection coating
Sensor	
Sensor type	CMOS
Resolution	5 Megapixel (2592x1944)
Maximum frame rate	30 fps
Compatibility	
Interface	USB 2.0
Operating system	Windows 7, 8, 10 & 11, MacOS 10.9 and up

Software	DinoCapture 2.0 (Windows), DinoXcope (Mac OS)
Supported image formats (Windows)	BMP, GIF, PNG, JPG, TIF, RAS, PNM, TGA, PCX, MNG, WBMP, JP2, JPC, PGX
Supported video formats (Windows)	WMV, FLV, SWF
Supported image formats (MacOS)	JPEG, PNG
Supported video formats (MacOS)	MOV (max 1.3MP)
Imaging standards	DirectShow, UVC
Wifi	Wireless-ready, requires the WF-10 WiFi streamer (optional)
Housing	
Housing material	Metal housing
Dimensions	10.5cm (L) x 3.2cm (D)
Weight	137g
Cable length	1.8m
Features	
Special feature	Automatic Magnification Reading (AMR), Extended Dynamic Range (EDR), Extended Depth of Field (EDOF), Flexible LED Control (FLC)
Measurement	Yes
Calibration	Yes
Microtouch sensor	Yes
ESD safe	Yes
Information	
Package contents	Microscope, carry pouch, software CD, calibration target, user manual, N3C-O- Open cap, N3C-C- Closed cap, N3C-D- Diffuser cap, N3C-E- Extension cap, N3C-L- Long cap, N3C-S- Side light cap
Warranty information	2 years European warranty
Regulatory approval	CE, FCC, ROHS
Price range	€1100,00 - €1300,00

Product Gallery