

AM4115T-DFRW

Brand:

Product Code: AM4115T-DFRW



Short Description

1.3 Megapixel resolution

620nm excitation and 650nm emission filter

20-220x magnification

To observe Cy5

Description

The Dino-Lite AM4115T-DFRW digital microscope is optimized for research and viewing fluorescent objects by using 620nm excitation LEDs. It has a long-pass type 650nm emission filter that is designed to observe a wide range of fluorophores with applications in developmental biology, pathology, and anatomy.



USB 2.0



Measurement
functionality



Adjustable
~20x-220x



Magnification
lock



1.3 megapixel
1280 x 1024



Fluorescent
LEDs ~ 620nm



Standard working
distance

The AM4115T-DFRW has the capability of switching the light source from the 620nm


excitation LEDs to the white LED, which is convenient for locating the object and obtaining an easy focus. All of these factors make the AM4115T-DFRW the smallest fluorescence microscope in the world, capable of many applications and visualizing a vast range of fluorescence.

» *Special models available on request.*

Working distance/field of view/depth of field

MAGNIFICATION RATE	working distance	FIELD OF VIEW(X)	FIELD OF VIEW(Y)	DEPTH OF FIELD
Listed values may differ slightly	* Without front cap			Unit = mm
20	60.2	19.5	15.6	2.5
30	33.5	13.0	10.4	1.8
40	20.9	9.8	7.8	1.5
50	13.9	7.8	6.3	-
60	9.7	6.5	5.2	-
70	7.1	5.6	4.5	1.0
80	5.5	4.9	3.9	-
90	4.5	4.3	3.5	-
100	4.1	3.9	3.1	-
110	4.0	3.6	2.8	-
120	4.1	3.3	2.6	-
130	4.5	3.0	2.4	-
140	5.0	2.8	2.2	-
150	5.6	2.6	2.1	-
160	6.3	2.4	2.0	-
170	7.1	2.3	1.8	-
180	8.0	2.2	1.7	-
190	8.9	2.1	1.6	-
200	9.9	2.0	1.6	-
210	10.9	1.9	1.5	-
220	11.9	1.8	1.4	0.1

Special Light

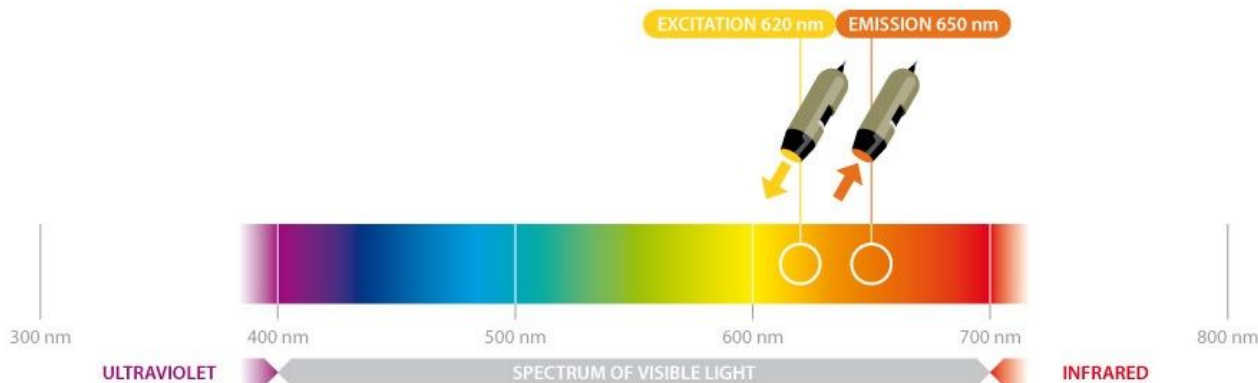

Specification

Lighting	
Light/ LED type	620nm excitation + white
Number of LEDs	7 FL + 1 White
LED on/off switchable:	Yes
Infrared filter	No
Diffuser available	No
Emission filter	650nm, LP
Polarizer	No
Fluorophores	Cy5
Optics	
Magnification	20x ~ 220x
Macro zoom	No
Working distance	Standard
Lens type	Glass with anti-reflection coating
Sensor	
Sensor type	CMOS
Resolution	1.3 Megapixel (1280x1024)
Maximum frame rate	30 fps
Compatibility	
Interface	USB 2.0
Operating system	Windows 7, 8, 10 & 11, MacOS 10.9 and up
Software	DinoCapture 2.0 (Windows), DinoXcope (Mac OS)
Supported image formats (Windows)	BMP, GIF, PNG, JPG, TIF, RAS, PNM, TGA, PCX, MNG, WBMP, JP2, JPC, PGX
Supported video formats (Windows)	WMV, FLV, SWF
Supported image formats (MacOS)	JPEG, PNG
Supported video formats	MOV

(MacOS)	
Imaging standards	DirectShow, UVC
Wifi	Wireless-ready, requires the WF-10 WiFi streamer (optional)
Housing	
Housing material	Composite/ plastic housing
Magnification lock	Yes
Dimensions	10.5cm (L) x 3.2cm (D)
Weight	100g
Cable length	1.8m
Features	
Special feature	No
Measurement	Yes
Calibration	Yes
Microtouch sensor	Yes
ESD safe	No
Information	
Package contents	Microscope, carry pouch, software CD, calibration target, user manual
Warranty information	2 years European warranty
Regulatory approval	CE, FCC, ROHS
Price range	€700,00 - €900,00