

RK-10A

Brand: Dino-Lite
Product Code: RK-10A

Short Description

High precision stand with quick release

Description

The Dino-Lite RK-10A is a sturdy and stable high-end stand solution that can be used with all the professional Dino-Lite models. It is constructed of resilient stainless steel and lightweight aluminum and offers a very precise fine-focus adjustment as well as a quick release function.

The quick release button allows for fast and convenient vertical movement and the fine-focus adjustment knob is especially useful at high magnification. This model also includes an adjustable safety stop to prevent the microscope from touching the specimen. Its adjustable horizontal arm has a 15cm range with 360-degree rotation and the additional arm (RK-10-EX included) increases its effective range to observe larger objects. With these features, the RK-10A is the ideal accessory for precise orientation and positioning for maximum productivity and ease of use in a wide range of environments including industrial, scientific and laboratory settings.


[AVAILABLE EXTENSIONS FOR THE RK-10A STAND](#)

RK-10-EX



additional horizontal arm (included with RK-10A)

RK-10-FX

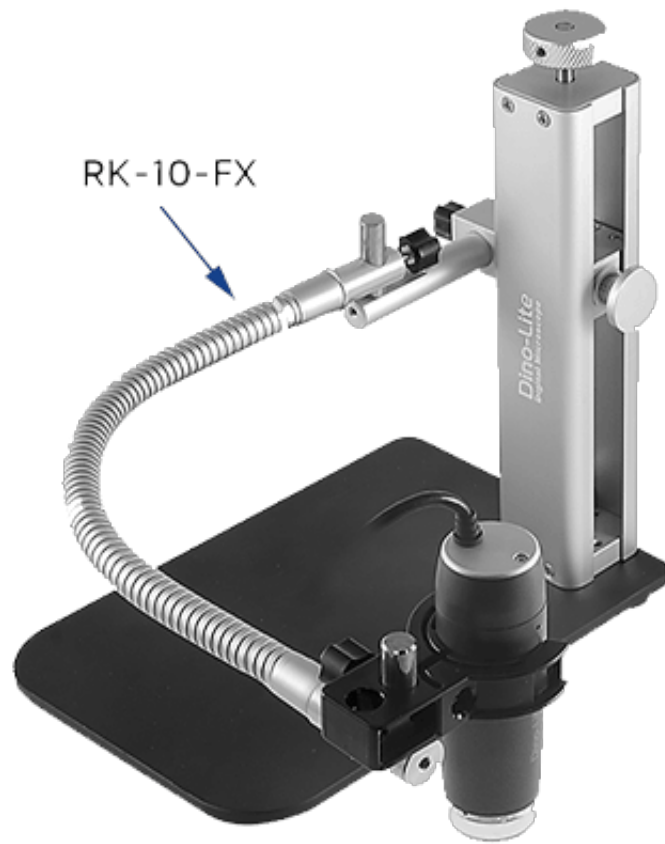

flexible arm extension

RK-10-VX



vertical arm extension

RK-10-PX



(Temporarily unavailable)

Specification

Compatibility	
Compatible with:	All Dino-Lite models
Housing	
Dimensions	15cm (W) x 22cm (L) x 23cm (H)
Weight	1.10 kg
Features	
ESD safe	Yes

Information

Package contents	Stand, HD-P1 holder, RK-10-EX horizontal arm
Warranty information	2 years European warranty
Accessory range:	Professional stands
Price range	€300,00 - €450,00

Product Gallery



